# **Symphony and Synapse LabVIEW VI User Manual**

# Table of Contents

Table of Contents .....	2
Overview: .....	3
Minimum Requirements: .....	3
User_Lib.LLB .....	4
Spectral Mode.vi .....	4
Image Mode.vi .....	5
Version number.vi .....	6
Glob_Lib.LLB .....	7
Read Detector uniq-id.vi .....	7
Initialize CCD.vi .....	8
Init ADC & Gain.vi .....	10
Setup CCD.vi .....	11
Spectral Reformat CCD area.vi .....	14
Image Reformat CCD area.vi .....	15
Acquire image.vi .....	16
Acquire spectral.vi .....	18
Acquire.vi .....	20
Read Temperature.vi .....	21
MultiAccu – Acquire Spectral.vi .....	23
MultiAccu – Configure.vi .....	25
MultiAccu - Disable.vi .....	27
MultiAccu - Get Operating Mode.vi .....	29
MultiAccu - Is Supported.vi .....	31
Trigger - Enable.vi .....	33
Trigger - Get Supported Input.vi .....	35
Trigger - Disable.vi .....	37
warning.vi .....	39
Appendix A:     Getting Started With HJY's LabVIEW VIs .....	40
Configuring an Instrument—REQUIRES SDK or SYNERJY .....	40
Configuring a Mono .....	40
Configuring a Single Channel Detector .....	40
Configuring a Multi-Channel Detector .....	41
Configuring a Filter Wheel .....	41
Initializing an Instrument .....	42
Appendix B:     Troubleshooting and Relinking HJY's LabVIEW VIs .....	44
Overview .....	44
Running on a Win 7 (32-bit and 64-bit) PC .....	47
Relinking VIs with ActiveX Controls: Examples .....	49
Service Policy .....	53
Return Authorization .....	54

---

## Overview:

This set of LabVIEW VIs is designed to run a Symphony, Synapse or Syncerity CCD system and provide the sub-VIs necessary to incorporate the Symphony or Synapse CCD into a full LabVIEW application.

The LabVIEW library, USER\_LIB.LLB contains two upper level VIs to control a CCD and collect data. These are meant as examples for a user to follow and modify for using a Symphony CCD in their application.

GLOB\_LIB.LLB contains the lower level sub-VIs that performs individual functions with the CCD. The individual VIs are used for the examples in USER\_LIB.LLB and documented in the manual.

## Minimum Requirements:

Windows 2000, XP or Windows 7 (See Appendix B for instructions on running with Windows7)

LabVIEW 7 or Higher

Installed Symphony CCD system

Installed Symphony COM component – This can be installed either using SynerJY software or the SynerJY SDK.


## User\_Lib.LLB

These are top level VIs that can be used as examples and a foundation for more advanced applications.

### *Spectral Mode.vi*

This is an example program to show how to use the sub-VIs in an application. This can be used to collect a spectrum or a series of spectra from different CCD Areas.


#### Front Panel
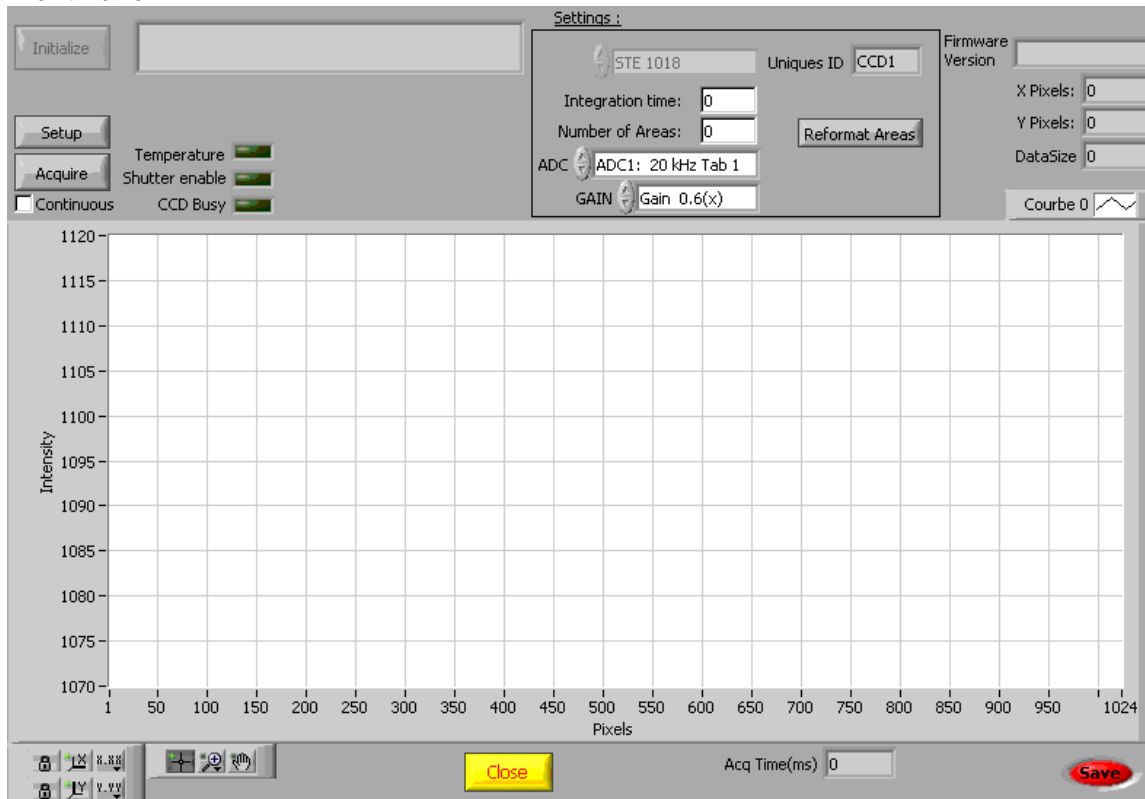
## Image Mode.vi

This is an example program to show how to use the sub-VIs in an application. This can be used to collect a spectrum or a series of spectra from different CCD Areas.

### Front Panel

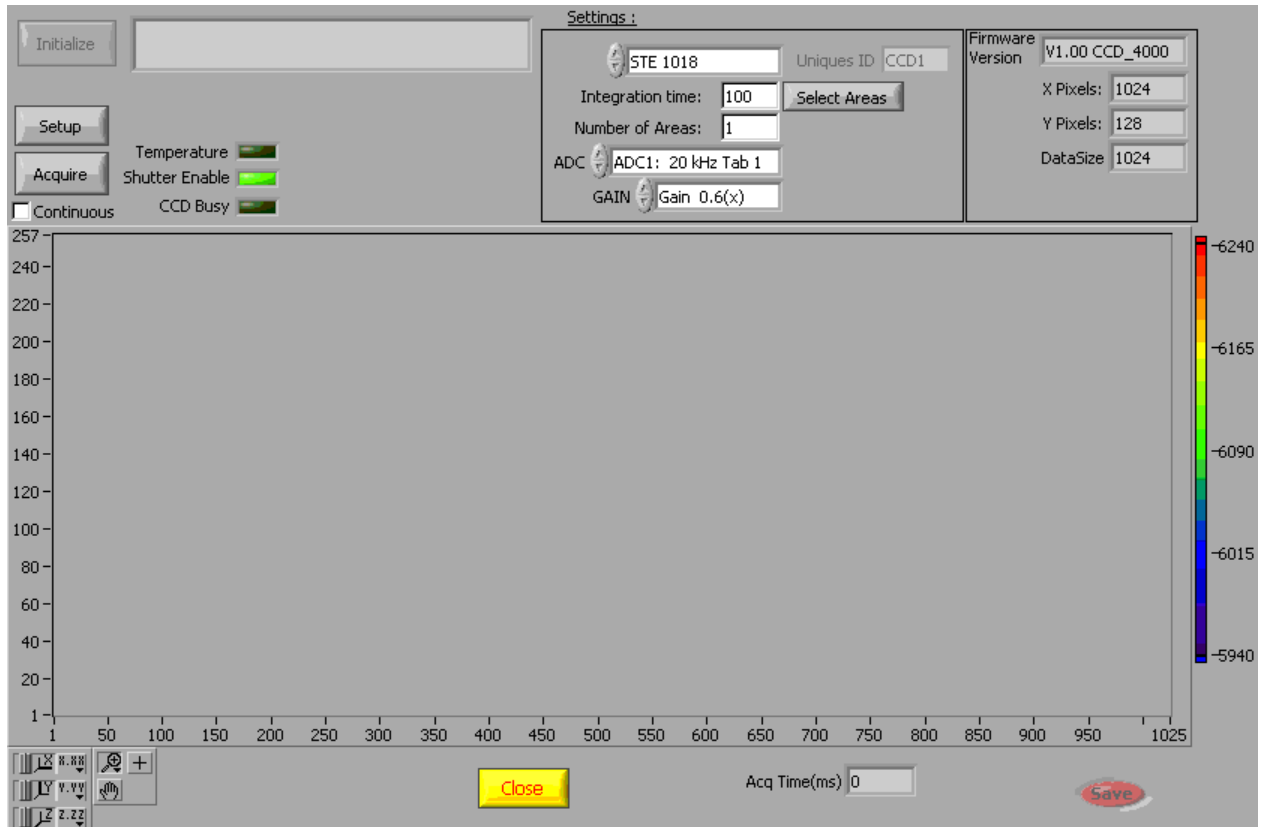---

## ***Version number.vi***

This is a Top Level VI that holds Version information about this revision of the Symphony LabVIEW VIs.

### **Connector Pane**

### **Front Panel**

---


## Glob\_Lib.LLB

These are the low level sub-VIs that perform different functions with the Symphony CCD

### *Read Detector uniq-id.vi*

This sub-VI will read all of the CCD information from the computer's registry. Registry settings are put into the registry on Hardware Configuration during installation of SynerJY software or the LabVIEW VIs.
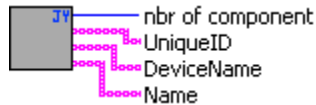
#### Connector Pane

#### Front Panel


#### Controls and Indicators

-  **nbr of component** : [Output] Total number of Detector UniqueIDs found in the registry
-  **UniqueID** [Output] Array of UniqueIDs
-  **DeviceName** : [Output] Array of Device Names corresponding to the same index in Unique ID Array
-  **Name** : [Output] Array of Specific Names corresponding to the same index in Unique ID Array


## Initialize CCD.vi

Sub VI to Initialize the CCD

### Connector Pane
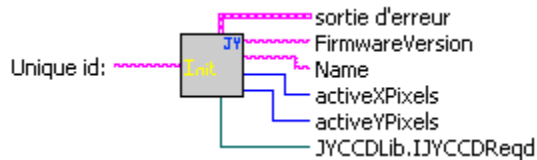
### Front Panel
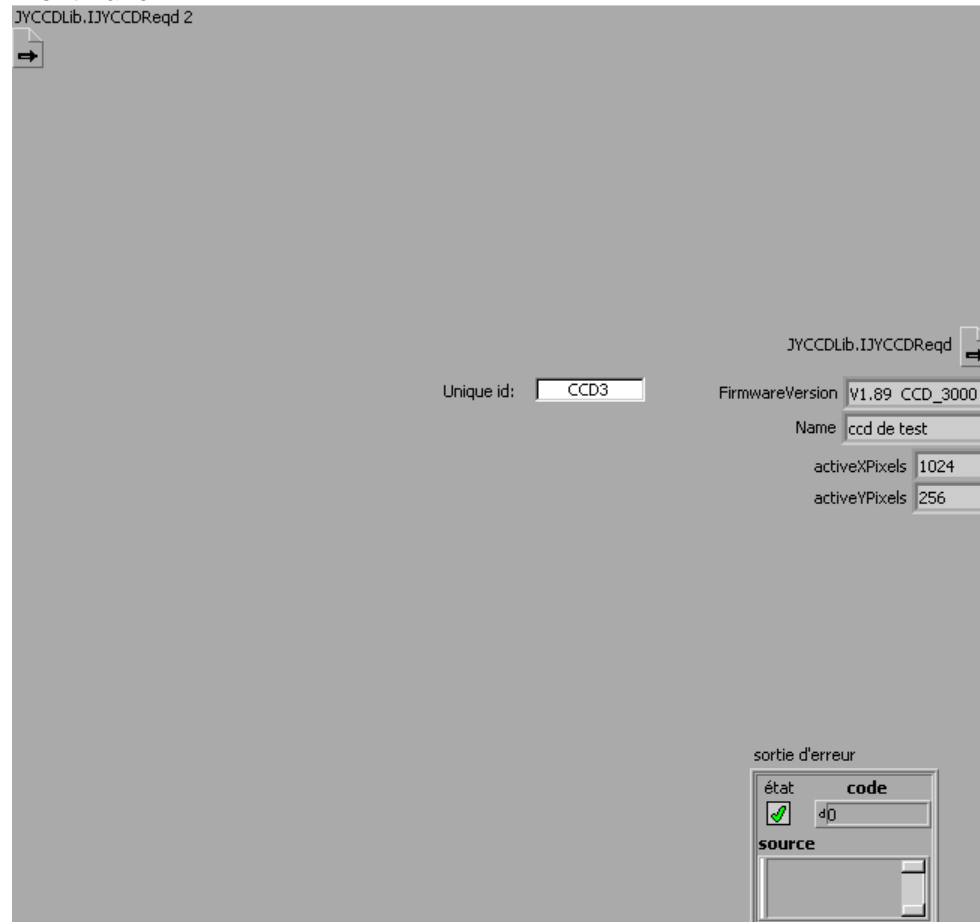
### Controls and Indicators

- Unique id:** [Input] - This is the Unique ID of the CCD. This information can be found by using the Read Detector uniq-id.vi VI
- activeYPixels:** [Output] Total number of Y pixels of the CCD
- ActiveXPixels:** [Output] Total number of X pixels of the CCD
- Name:** [Output] Device Name of CCD
- Firmware Version:** [Output] Firmware version of the Symphony CCD Controller
- Output error (Sortie d'erreur)** The **error out** cluster can pass error information to other VIs in the software. Use this information to decide if any functionality should be

---

bypassed in the event of errors from other VIs.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

**status** The **status** boolean is either TRUE (X) for an error, or FALSE (checkmark) for no error or a warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

**code** The **code** input identifies the error or warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

**source** The **source** string describes the origin of the error or warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.


**JYCCDLib.IJYCCDReqd**


Link to the CCD COM component

## Init ADC & Gain.vi

### Connector Pane
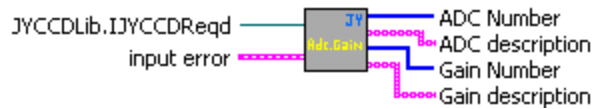
### Front Panel
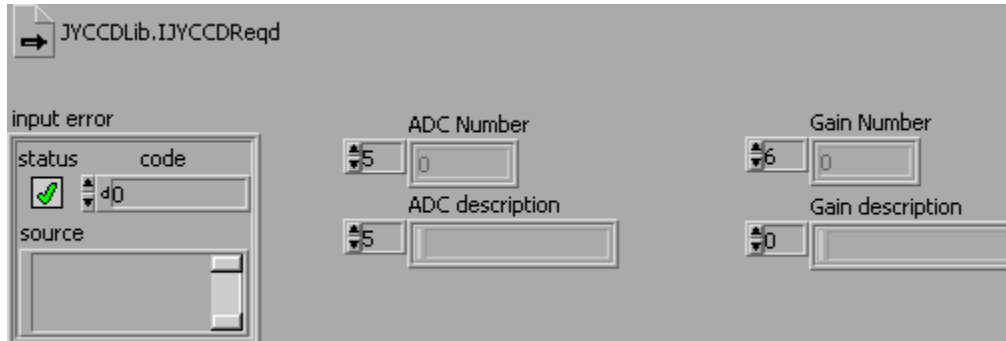
### Controls and Indicators

**JYCCDLib.IJYCCDReqd**

**Input error (Entree d'erreur)** The **error in** cluster can accept error information wired from VIs previously called. Use this information to decide if any functionality should be bypassed in the event of errors from other VIs.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

**status** The **status** boolean is either TRUE (X) for an error, or FALSE (checkmark) for no error or a warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

**code** The **code** input identifies the error or warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

**source** The **source** string describes the origin of the error or warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

**ADC description** – [Output] Array of ADC Descriptions from the CCD Controller

**ADC Number** – [Output] Array of ADC Values from the CCD Controller


**Gain Number** – [Output] Array of Gain Descriptions from the CCD Controller

## Setup CCD.vi

VI to setup a CCD for an acquisition

### Connector Pane
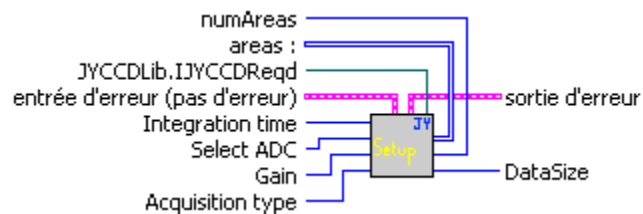
### Front Panel

The front panel of Setup CCD.vi includes the following controls and indicators:

- JYCCDLib.IJYCCDReqd**: A button in the top left corner.
- DataSize**: A numeric control set to 0.
- Integration time**: A numeric control set to 10.
- Select ADC**: A dropdown menu set to JY\_ADC\_16\_BIT.
- Gain**: A dropdown menu set to gain 1.
- Acquisition type**: A dropdown menu set to JYMCD\_ACQ\_FORMAT\_SCAN.
- numAreas**: A numeric control set to 2.
- areas:** A table with 2 rows and 5 columns. The first column contains values 1 and 1. The second column contains values 1024 and 256. The other three columns contain 0.
- entrée d'erreur (pas d'erreur)**: An indicator with a green checkmark and a code field set to 0.
- sortie d'erreur**: An indicator with a green checkmark and a code field set to 0.

### Controls and Indicators

**Gain** : [Input] Set the Gain Level for the CCD Acquisition

**Integration time** : [Input] Set the Integration Time for the Acquisition

---

 **Select ADC** : [Input] Select the ADC for the acquisition

 **numAreas** : [Input]

 **Acquisition type** [Input]

 **JYCCDLib.IJYCCDReqd**

This is a link to the CCD COM object functions

 **Input error (Entree d'erreur)** The **error in** cluster can accept error information wired from VIs previously called. Use this information to decide if any functionality should be bypassed in the event of errors from other VIs.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

 **status** The **status** boolean is either TRUE (X) for an error, or FALSE (checkmark) for no error or a warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

 **code** The **code** input identifies the error or warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

 **source** The **source** string describes the origin of the error or warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

 **areas** :

 **Output error (Sortie d'erreur)** The **error out** cluster can pass error information to other VIs in the software. Use this information to decide if any functionality should be bypassed in the event of errors from other VIs.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

 **status** The **status** boolean is either TRUE (X) for an error, or FALSE (checkmark) for no error or a warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

 **code** The **code** input identifies the error or warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

 **source** The **source** string describes the origin of the error or warning.


The pop-up option **Explain Error** (or Explain Warning) gives more information

---

about the error displayed.

## Spectral Reformat CCD area.vi

### Connector Pane
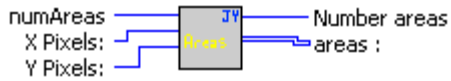
### Front Panel

areas			
0	1	0	0
0	1	0	0
1	0	0	0
1024	0	0	0
256	0	0	0
1	0	0	0
256	0	0	0

### Controls and Indicators:

**numAreas** [Input]

Total Number of Areas defined

**X Pixels** [Input]

Total Number of pixels in the X direction of the CCD

**Y Pixels** [Input]

Total Number of pixels in the Y direction of the CCD

**Areas** [Output]

**Area Number:** specific number of area for this set of parameters

**X Start:** First pixel in X dimension

**Y Start:** First pixel in y dimension

**X Size :** Total number of pixels in X direction

**Y Size :** Total number of pixels in Y direction


**X binning :** Number of pixels to sum together to get 1 data value in X direction

**Y binning:** Number of pixels to sum together to get 1 data value in Y direction.  
For spectral information, this value needs to be equal to the Y size.

## Image Reformat CCD area.vi

Sub-VI to define an area for acquiring an image from the CCD

### Connector Pane
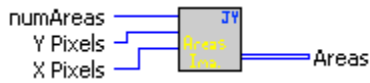
### Front Panel

Area number	X Start	Y Start	X Size	Y Size	X Bining	Y Bining
1	0	0	1024	256	10	10

### Controls and Indicators

**numAreas** [Input]

Total Number of Areas defined

**X Pixels** [Input]

Total Number of pixels in the X direction of the CCD

**Y Pixels** [Input]

Total Number of pixels in the Y direction of the CCD

**Areas**

**Area Number:** specific number of area for this set of parameters

**X Start:** First pixel in X dimension

**Y Start:** First pixel in y dimension

**X Size :** Total number of pixels in X direction

**Y Size :** Total number of pixels in Y direction


**X binning :** Number of pixels to sum together to get 1 data value in X direction

**Y binning :** Number of pixels to sum together to get 1 data value in Y direction


## Acquire image.vi

This Sub-VI starts the acquisition of data from the Symphony CCD and returns a 3D array with X Data Value (Pixel), Y Data Value (Pixel), and Intensity

### Connector Pane
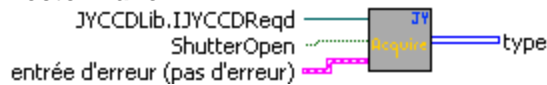
### Front Panel
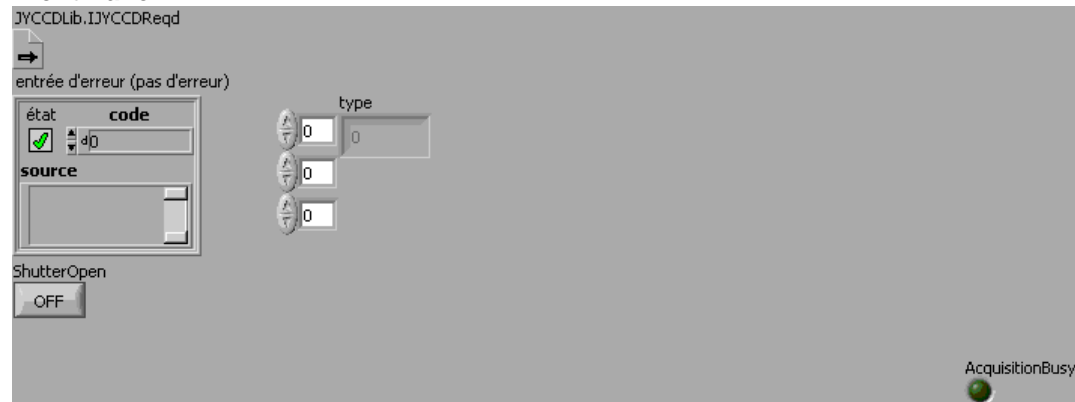
### Controls and Indicators

#### JYCCDLib.IJYCCDReqd

This is a link to the CCD COM object functions

 **Input error (Entree d'erreur)** The **error in** cluster can accept error information wired from VIs previously called. Use this information to decide if any functionality should be bypassed in the event of errors from other VIs.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

 **status** The **status** boolean is either TRUE (X) for an error, or FALSE (checkmark) for no error or a warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

 **code** The **code** input identifies the error or warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

 **source** The **source** string describes the origin of the error or warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

#### ShutterOpen

True to open the shutter during acquisition

---

**[I32]**


False to keep the shutter closed during an acquisition  
**type**

data information from the image being acquired.


## Acquire spectral.vi

This Sub-VI starts the acquisition of data from the Symphony CCD and returns a 2D array with X Data Value (Pixel) and Intensity

### Connector Pane
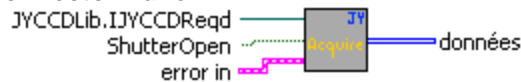
### Front Panel
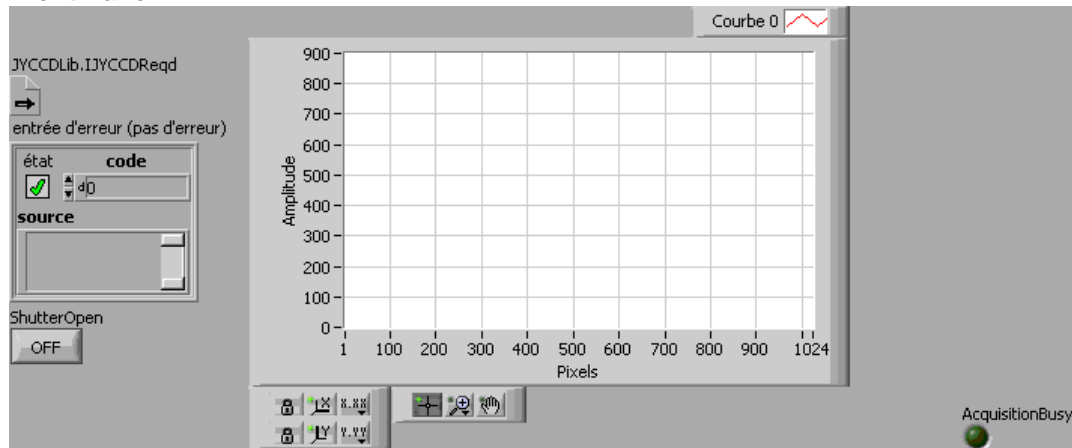
### Controls and Indicators

**JYCCDLib.IJYCCDReqd**

This is a link to the CCD COM object functions

**Input error (Entree d'erreur)** The **error in** cluster can accept error information wired from VIs previously called. Use this information to decide if any functionality should be bypassed in the event of errors from other VIs.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

**status** The **status** boolean is either TRUE (X) for an error, or FALSE (checkmark) for no error or a warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

**code** The **code** input identifies the error or warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

**source** The **source** string describes the origin of the error or warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

**ShutterOpen**

---

**I32**


True to open the shutter during acquisition  
False to keep the shutter closed during an acquisition  
**type**

data information from the image being acquired.


## Acquire.vi

This Sub-VI starts the acquisition of data from the Symphony CCD and returns a 3D array with X Data Value (Pixel), Y Data Value (Pixel), and Intensity. This VI has the same function and Acquire Image.

### Connector Pane
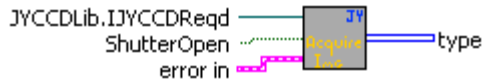
### Front Panel
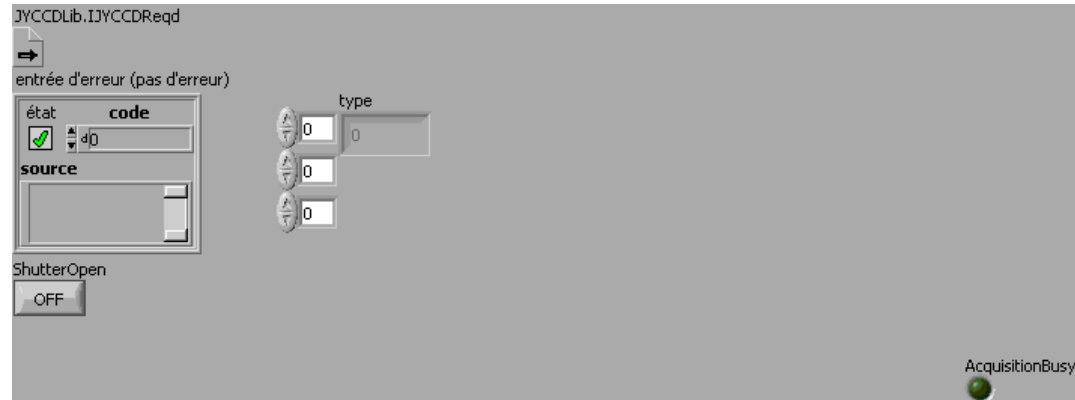
### Controls and Indicators

#### JYCCDLib.IJYCCDReqd

This is a link to the CCD COM object functions

 **Input error (Entree d'erreur)** The **error in** cluster can accept error information wired from VIs previously called. Use this information to decide if any functionality should be bypassed in the event of errors from other VIs.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

 **status** The **status** boolean is either TRUE (X) for an error, or FALSE (checkmark) for no error or a warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

 **code** The **code** input identifies the error or warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

 **source** The **source** string describes the origin of the error or warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

 **ShutterOpen** : Input – Set to True to open the shutter during acquisition. Set to False to keep shutter closed during acquisition.

 **Type**


---

Array of data reported back from detector. Data is in the format of X, Y, Intensity


## ***Read Temperature.vi***

Sub-VI to read the CCD's current temperature

### **Connector Pane**
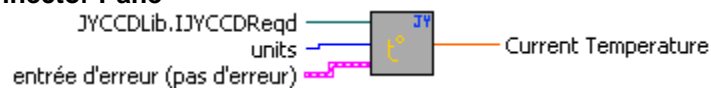
### **Front Panel**
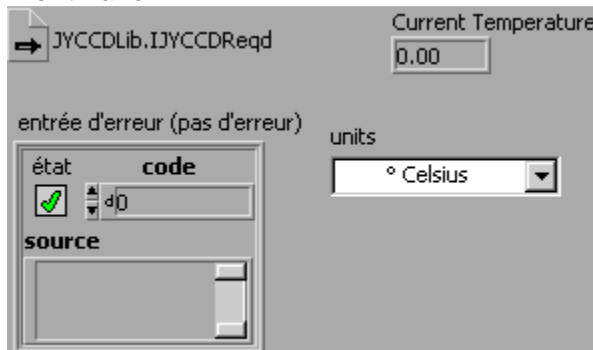
### **Controls and Indicators**

**JYCCDLib.IJYCCDReqd**

This is a link to the CCD COM object functions

**Input error (Entree d'erreur)** The **error in** cluster can accept error information wired from VIs previously called. Use this information to decide if any functionality should be bypassed in the event of errors from other VIs.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

**status** The **status** boolean is either TRUE (X) for an error, or FALSE (checkmark) for no error or a warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

**code** The **code** input identifies the error or warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

**source** The **source** string describes the origin of the error or warning.

The pop-up option **Explain Error** (or Explain Warning) gives more information about the error displayed.

**units** [Input]: Select Celsius, Fahrenheit, or Kelvin temperature units


**Current Temperature** [Output]: Shows the current CCD Temperature

---


## MultiAccu – Acquire Spectral.vi

Sub-VI to acquire data in Multi-Accumulation Mode.

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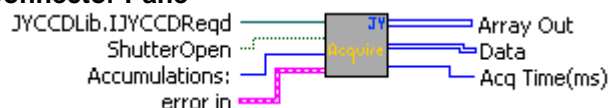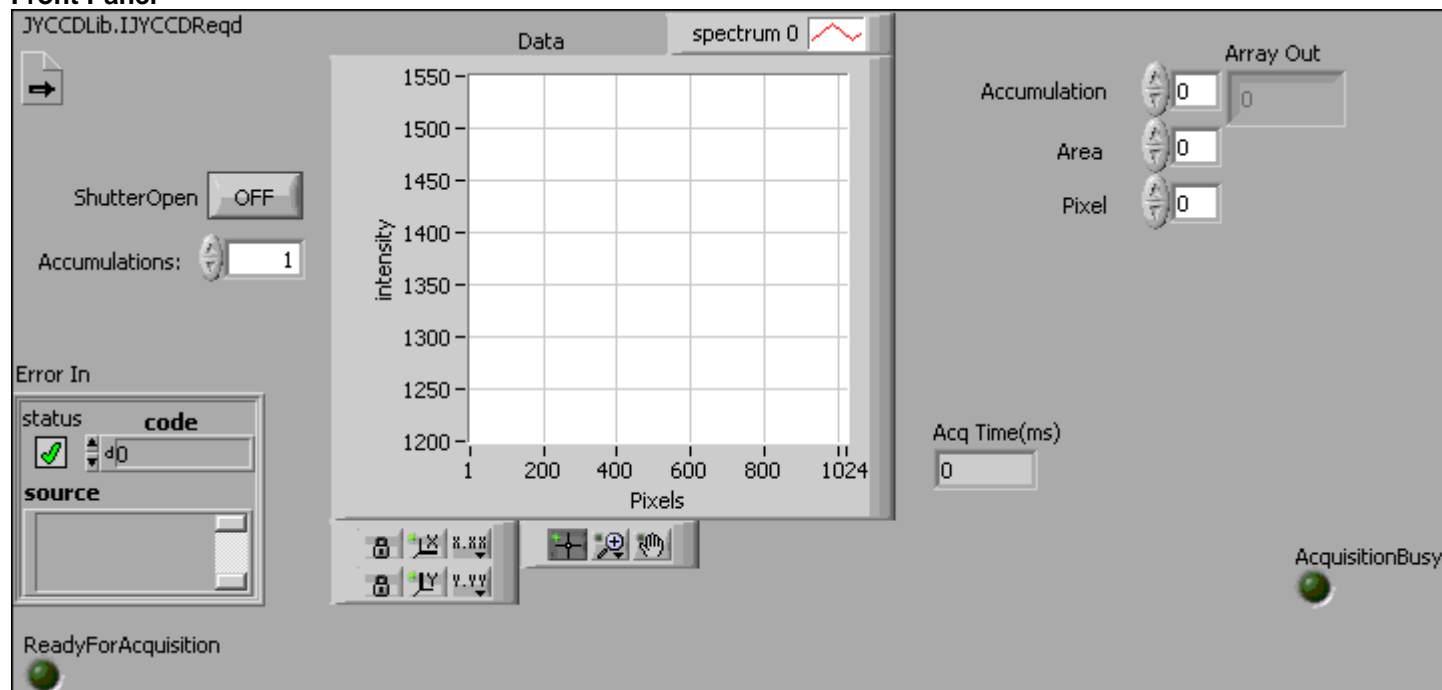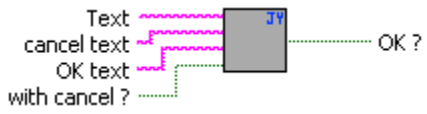
OK ?

---

## Appendix A:      Getting Started With HJY's LabVIEW VIs

---

HORIBA Scientific provides full sets of dedicated LabVIEW VIs for controlling its Monos, Detectors, and Filter Wheels. The modular design of the VIs makes it easy for them to interface with HORIBA Jobin Yvon instruments and with each other. The HORIBA Scientific LabVIEW VIs require LabVIEW Version 7 or higher.

In general, it is recommended that instruments be configured using SynerJY, SynerJY-SDK, USB Spectrometer Control, or in certain cases Lynear (for certain CCDs). This document describes a suggested procedure for establishing communications with an instrument.

### Configuring an Instrument—REQUIRES SDK or SYNERJY

A valid registry configuration is required to control the instrument. The **JYAddDevice.vi** provides a simple way to configure a device. Open and run the VI.

#### Configuring a Mono

1. Select **Monochromator** from the list and click **Next**.
2. Select a **Monochromator type** from the list and click **Next**.
3. If there is more than one type of communications available for the mono, you will be able to choose. If **GPIB** is selected make sure to enter the valid GPIB address in the **Port Number** field. If **Serial** is selected make sure to enter the correct **COM Port Number** and the correct values for the **Baud Rate**, **Data Bits**, **Stop Bits** and **Parity** bits. Click **Next**.
4. Assign any display name you like in the field provided. Click **Next**.
5. Enter the correct values for **Groove Density**, **Blaze**, **Mirrors** and **Shutter**. Click **Next**.
6. Click **Apply**.

#### Configuring a Single Channel Detector

1. Select **Detector** from the list and click **Next**.
2. Select **Single Channel Detector** from the list and click **Next**.
3. Select the **Detector Type** (Lock-in, SpectrAcq2) from the **Special Types** combo box.
4. Make the correct selections for **Options** and **Communications Parameters** for the selected type. Click **Next**.
5. Assign any display name you like in the field provided. Click **Next**.

- 
6. Click **Apply**.

**Note:** Although the SCD VI toolkit VIs will allow simple data acquisition for the SR810 and similar lock-in amplifiers, it is recommended that the user download the manufacturers VIs.


## Configuring a Multi-Channel Detector

1. Select **Detector** from the list and click **Next**.
2. Select **Multi Channel Detector** and click **Next**.
3. Select the **Controller Type** from the **Special Types** combo box. Make the correct selections from the choices provided. Click **Next**.
4. Assign any display name you like in the field provided. Click **Next**.
5. If a Symphony was selected you will be required to enter the path to the table files. If you have selected a detector type that does not require table files, this screen will not be displayed. Click **Next**.
6. Click **Apply**.

## Configuring a Filter Wheel

1. Select **Accessory** from the list and click **Next**.
2. Select **Filter Wheel** and click **Next**.
3. Select the **Type** of Filter Wheel(Triax, USB) and set the appropriate **Communications Parameters**. Click **Next**.
4. Assign any display name you like in the field provided. Click **Next**.
5. Set the **Number of Positions** in the filter wheel and set the **Wavelength Info** for each filter wheel position. Click **Next**.
6. Click **Apply**.

The following page lists the VIs available for HORIBA Scientific filter wheels.
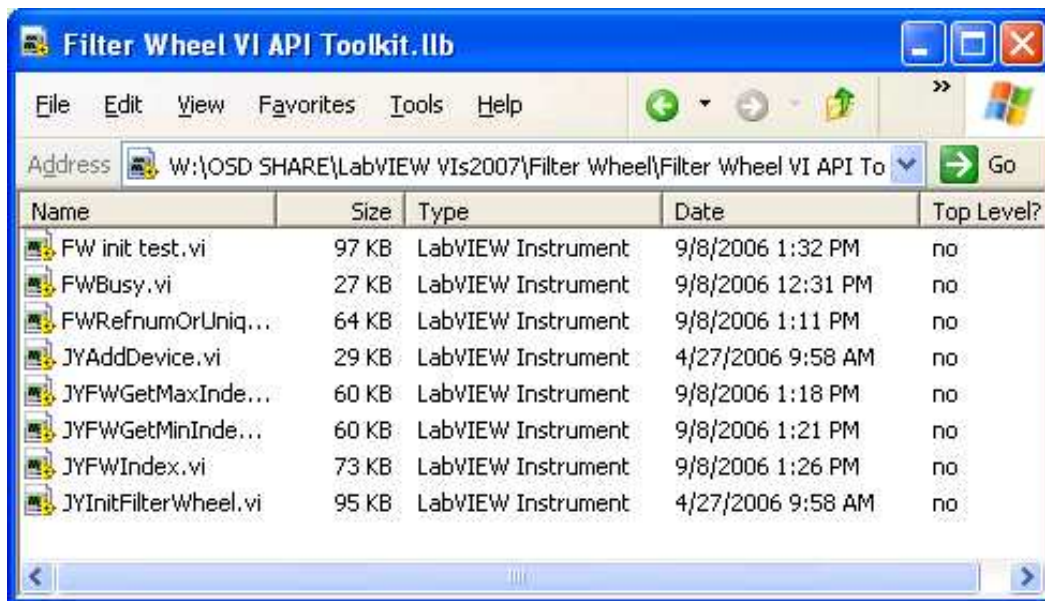
## Initializing an Instrument

This section describes how to initialize/establish communications with a monochromator. The procedure is identical for filter wheels and detectors.

A valid **UniqueID** is required to initialize a mono. A valid UniqueID is created for each configured device.

### *Configuring an Instrument before continuing with this section:*

To find the UniqueID for a configured mono use the GetMonosFromConfigBrowser.vi. This will provide two lists, a list of UniqueIDs and a list of Display Names. There is a one to one correspondence between the lists. The first UniqueID in the UniqueID list corresponds to the first Display Name in the Display Name List.

Open the **InitMono.vi**. Type the **UniqueID** in the space provided and run the VI. If you have provided the correct parameters when configuring the mono, the mono will go through its initialization sequence.

If the mono fails to initialize, check that a) the mono is powered on; b) the mono is connected to the PC; 3) the **Communications Parameters** (Serial, USB, Baud Rate, Port Number, etc) are correct.

Once you have successfully initialized the mono you are ready to start using the HORIBA Scientific LabVIEW drivers to create your own application.

---

**Note:** In order to control and acquire data with a Single Channel Detector, the SCD VI Toolkit VIs are required. This Toolkit is included with the HORIBA Jobin Yvon full installation version of the LabVIEW USB Monochromator VI. The VIs for multichannel detectors are separate and include a separate user manual. To control lock-in amplifiers, it is recommended to use the VIs provided by the lock-in amplifier manufacturer.

## Appendix B: Troubleshooting and Relinking HJY's LabVIEW VIs

---


HORIBA Scientific provides full sets of dedicated LabVIEW VIs for controlling its Monos, Detectors, and Filter Wheels. The modular design of the VIs makes it easy for them to interface with HORIBA Jobin Yvon instruments and with each other. The HORIBA Jobin Yvon LabVIEW VIs require LabVIEW Version 7 or higher.

### Overview


Before Starting the LabVIEW VIs for HORIBA Scientific CCDs, make sure SynerJY or the SynerJY-SDK has been installed on the computer (or Lynear for Sygnature detectors). For USB monos, USB Spectrometer Control will also work. It is easiest to create the mono or CCD through one of these HJY software packages.

If the VIs for Symphony, Synapse or Sygnature CCDs crash on the first execution of the system, the following procedure should be performed on each sub-VI in the HORIBA Jobin Yvon Symphony VI library, Global\_Lib.LLB. After all of the sub-VIs have been re-linked and saved, the Symphony VIs should work correctly.

1. Start LabVIEW and open the sub-VI
2. Click on the Window menu and select 'Show Block Diagram'
3. Locate the first ActiveX Property Node in the sub-VI

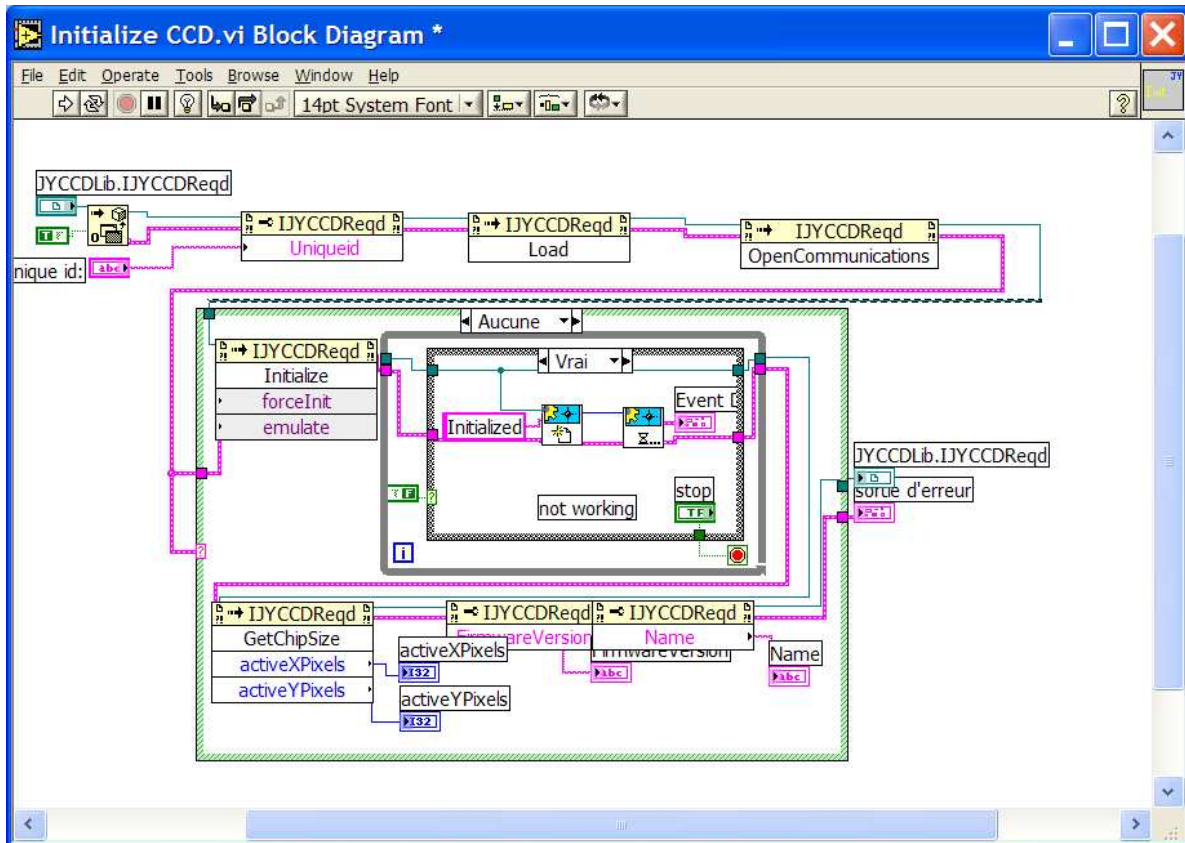4. Record the current property of the Property Node. Click on the Property Name to get the list of available properties and change the property to a different property than was originally listed.
5. Click on the Property Name again to get the list of available properties and change the property back to the original property that was listed.
6. Locate the first ActiveX Invoke Node in the sub-VI

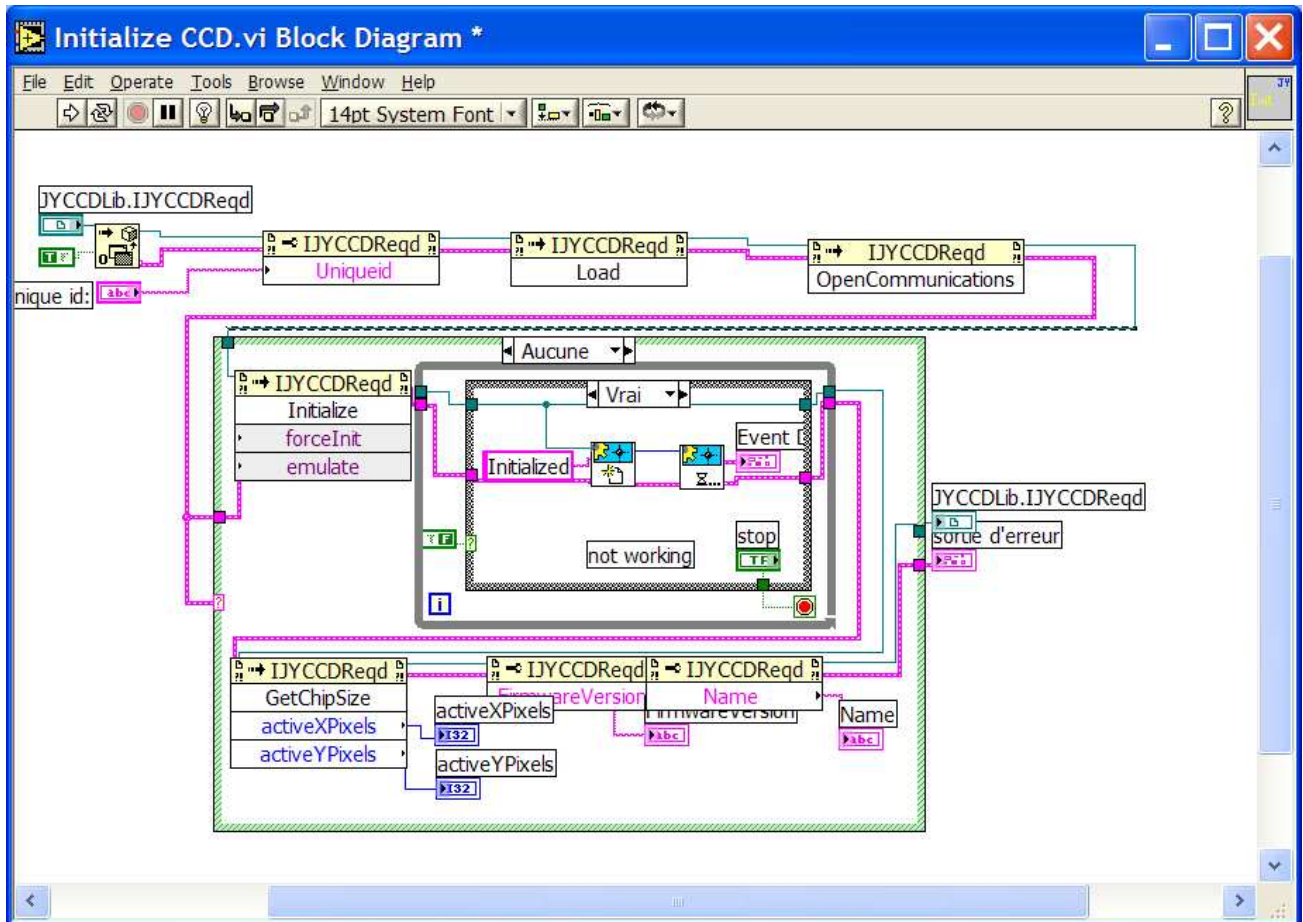7. Record the current method of the Invoke Node. Click on the Method Name to get the list of available properties and change the method to a different method than was originally listed.
8. Click on the Method Name again to get the list of available methods and change the method back to the original method that was listed.
9. Save the sub-VI.

If errors continue, turn on the Trace mode to determine where the error is occurring. If there is an Active X Invoke or Property Node or a sub-VI that has not been re-linked, perform the above procedure on that sub-VI.

Other items to check:

- The user has selected a valid Device and UniqueID in the top level VI.
- If the user has modified the VIs, make sure the user is passing the ActiveX Ref Number and Error cluster to each of the sub-VIs after the CCD has been initialized.

## Running on a Win 7 (32-bit and 64-bit) PC

If a crash/hang up occurs when attempting to communicate with a monochromator or other instrument using LabVIEW on a Windows 7 system running on a 32 bit or 64 bit platform, this most probably is due to something called Data Execution Prevention (DEP). Built into the Windows 7 framework, this prevents/monitors certain programs from using up too many resources which can be a symptom of a virus or other malware. To disable DEP, do the following:

- 1) Shut the computer off.
- 2) Start the computer pressing the F2 button (or Enter) to open the BIOS.
- 3) As this next part varies from manufacturer to manufacturer, please note the following:
  - a. On some computers, the heading of “Security” reveals a new heading of “Memory Protection”. Disable this feature to disable DEP.
  - b. On some computers, the heading of no-execute page-protection (NX) must be disabled (AMD processors)
  - c. On some computers, under Advanced→CPU→Execute Disable Bit must be disabled
  - d. The Execute BIT (XD) must be toggled to Disabled (Intel Processors).
- 4) Press F10 to save and quit restarting the computer. DEP should now be disabled and communications should not hang up with respect to the VI and instrument.

Please have good virus protection software running and up to date as this additional safety check is now disabled.

LabVIEW must be run in compatibility mode in Windows 7 to work correctly. To set compatibility:

- 1) Right-click the LabVIEW icon and choose Troubleshoot Compatibility.
- 2) You can either choose Windows XP Service Pack 2 manually or Windows can automatically decide compatibility.
- 3) Launch LabVIEW to be sure that the problem is fixed. If Windows asks if you want to allow LabVIEW to make changes to your computer, click yes.
- 4) Save the settings so that LabVIEW will be opened in compatibility mode in the future.

When using a detector with USB interface (such as Synapse, Symphony or Sincerity CCD and IGA array detectors), there are additional considerations. The primary difference between running USB detectors in Windows XP and Windows 7 is where in the registry the devices are logged. The locations will differ depending on operating system. The location can be verified through Windows regedit program. In Windows XP and Windows 7 32-bit, USB devices will appear in the registry at the following location:

HKEY\_LOCAL\_MACHINE\SOFTWARE\Jobin Yvon\Components\Detectors


In Windows 7 64-bit, USB detectors will appear in the registry at the following location:

HKEY\_LOCAL\_MACHINE\SOFTWARE\Wow6432Node\Jobin

Yvon\Components\Detectors


The difference in registry location becomes crucial in the sub-vi “read detector uniq-id.vi”

## Windows XP and Windows 7 32-bit

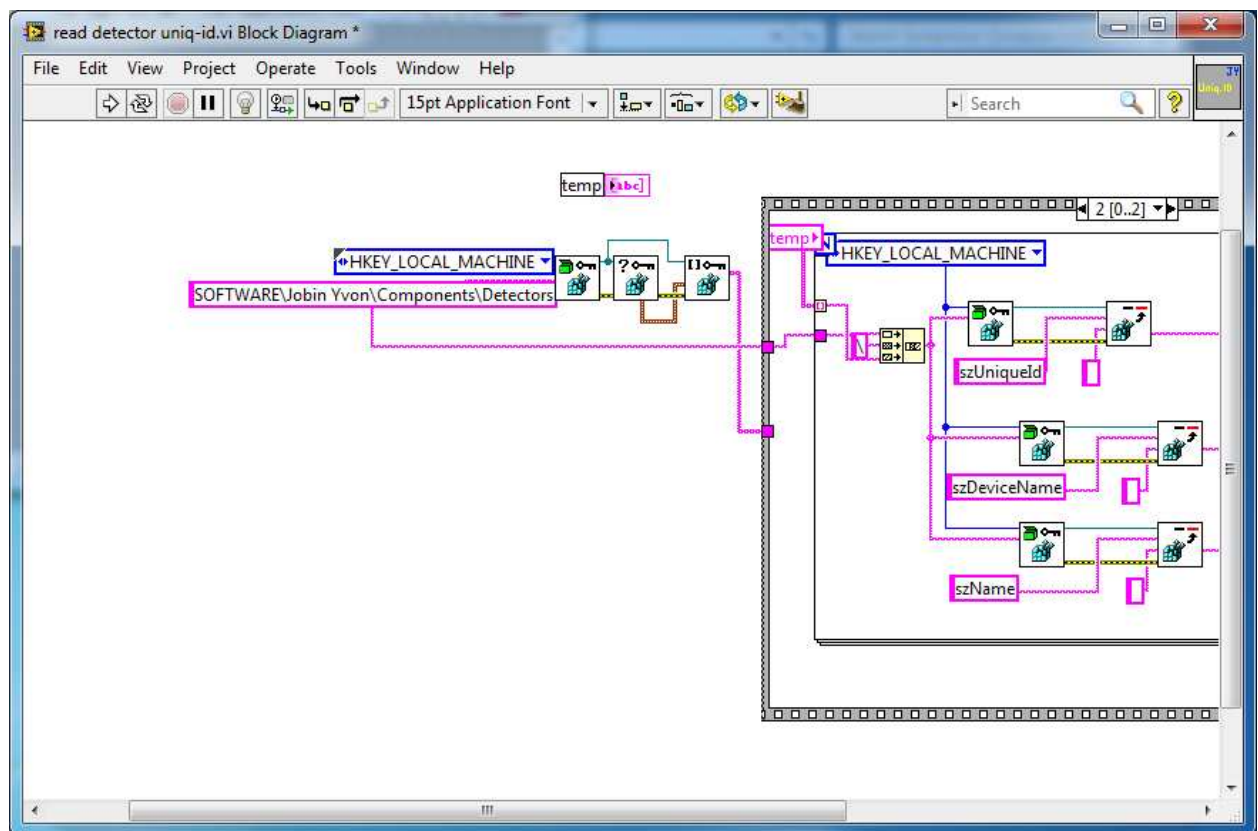
## Windows 7 64-bit
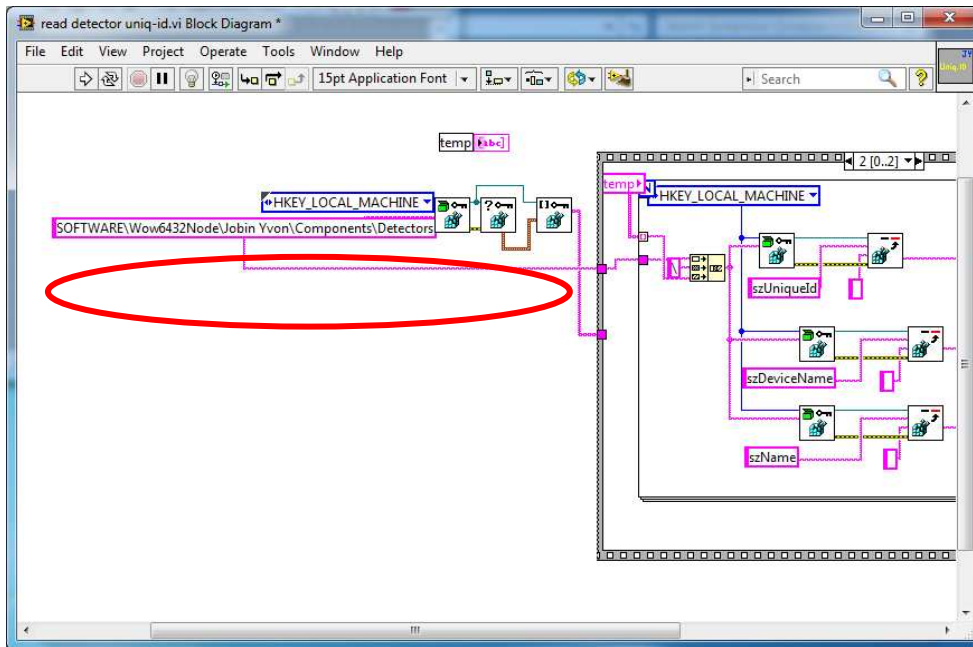
## Relinking VIs with ActiveX Controls: Examples

On occasion, a LabVIEW VI containing ActiveX controls may become broken or “unlinked” from the ActiveX control. Usually this happens if a newer version of the ActiveX control is installed on a PC or the VI is opened in a newer version of LabVIEW. In either case, the linkage between the LabVIEW Automation Ref Num and the ActiveX control must be re-established. The following procedure describes that process.


1. **OPEN THE VI**—Open the Front Panel of a “broken” VI. (A VI is broken when the arrow is cut in two parts.)

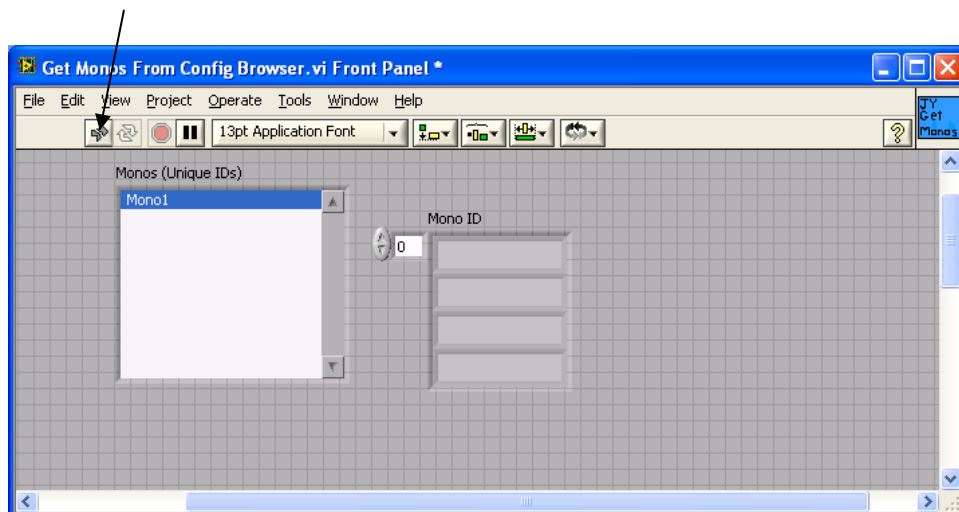
2. **LOCATE THE AUTOMATION REFNUM.** This symbol looks like a page with the upper right hand corner folded down. On the page is a black, right pointing arrow.

3. **Select The ActiveX Class**

Right click on the control and from the pop-up menu click on “Select ActiveX Class” and write down the ActiveX class that is selected. In the example above, it is JYMONOLib.IJYMonoReqd. (Note: all monochromator control functions will have the same class (JYInitMono, JYMirrorMove, JYMoveToWavelength, etc.) Now click on “Browse...” to bring up the “SelectObject From Type Library dialog”. Click the blue down arrow on the Type Library combo box to select the ActiveX class from the list. If the ActiveX class is not in the list add by clicking on the Browse.. button. After the ActiveX class is selected, one or more entries will appear in the lower box titled “Objects”. Select the desired ActiveX class and click OK. (In some versions of LabVIEW instead of a combo box there will be a list. First, select a different ActiveX class than the one you just wrote down, then reselect the original one and Save.). The first example shown is with a monochromator control VI (JYSlits.VI)

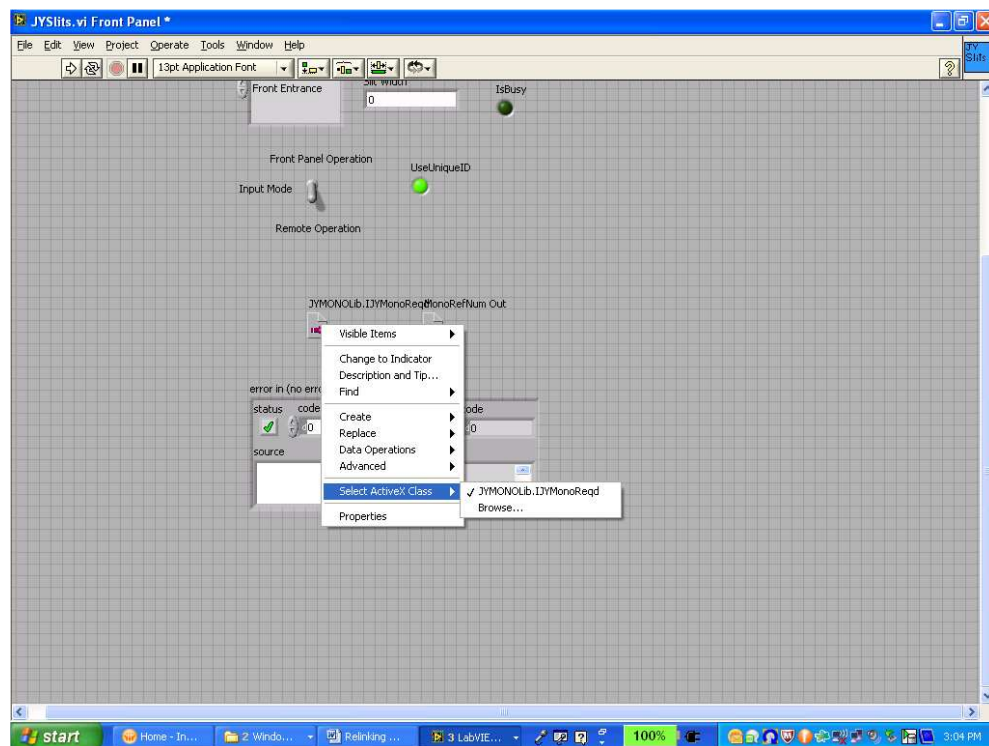
The second example is a config browser type VI using ActiveX class JYCONFIGBROWSERCOMPONENTLib.IJYConfigBrowserInterface
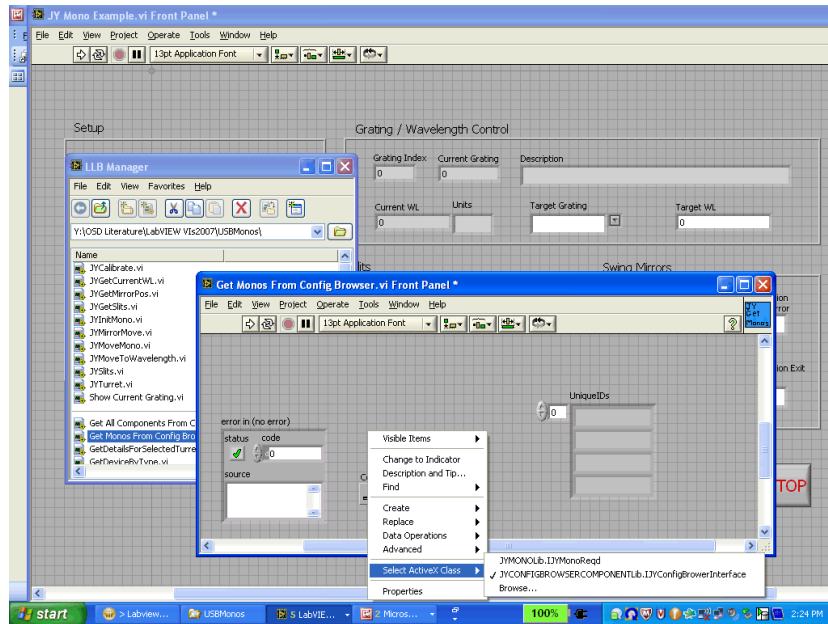
#### 4. Wire the Invoke and Property Nodes

Open the block diagram for the broken VI. All the invoke and property nodes wired to this Automation Refnum will now appear in light gray as if they are disabled. Place the mouse over the invoke node or property node. When the cursor changes to a pointing finger, left click to bring up the list of methods or properties associated with that node. Select the desired property or method from the pop up list. Repeat this procedure for each “broken” invoke node or property node on the block diagram. Save the VI.

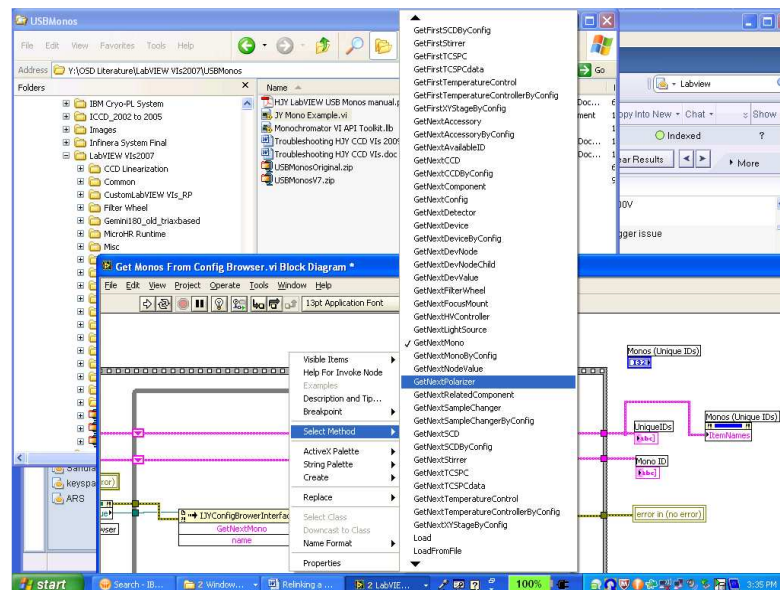- 5. Repeat this procedure for each VI that needs relinking.**

## Service Policy

---

If you need assistance in resolving a problem with your instrument, contact our Customer Service Department directly, or if outside the United States, through our representative or affiliate covering your location.

Often it is possible to correct, reduce, or localize the problem through discussion with our Customer Service Engineers.

All instruments are covered by warranty. The warranty statement is printed on the inside back cover of this manual. Service for out-of-warranty instruments is also available, for a fee. Contact HORIBA Jobin Yvon or your local representative for details and cost estimates.

If your problem relates to software, please verify your computer's operation by running any diagnostic routines that were provided with it. Please refer to the software documentation for troubleshooting procedures. If you must call for Technical Support, please be ready to provide the software serial number, as well as the software version and firmware version of any controller or interface options in your system. The software version can be determined by selecting the software name at the right end of the menu bar and clicking on "About." Also knowing the memory type and allocation, and other computer hardware configuration data from the PC's CMOS Setup utility may be useful.

In the United States, customers may contact the Customer Service department directly. From other locations worldwide, contact the representative or affiliate for your location.

### **In the USA:**

HORIBA Scientific  
3880 Park Avenue  
Edison, NJ 08820-3012  
Toll-Free: +1-866-Jobinyvon  
Tel: +1-732-494-8660 Ext. 160  
Fax: +1-732-549-5125  
Email:  
[service.jyus@horiba.com](mailto:service.jyus@horiba.com)  
[www.horiba.com/scientific](http://www.horiba.com/scientific)

### **In France:**

HORIBA Jobin Yvon  
S.A.S  
16-18 rue du Canal  
91165 Longjumeau cedex  
Tel: +33 (0) 1 64 54 13 00  
Fax: +33 (0) 1 69 09 07 21  
[www.jobinyvon.fr](http://www.jobinyvon.fr)

### **In Japan**

Horiba Ltd., JY Optical Sales Dept.  
Higashi-Kanda Daiji Building  
1-7-8 Higashi-Kanda, Chiyoda-ku  
Tokyo 101-0031  
Tel: +81 (0) 3 3861 8231  
[www.jyhoriba.jp](http://www.jyhoriba.jp)

### **Worldwide: 1-877-JYHoriba**

**China:** +86 (0) 10 6849 2216

**Germany:** +49 (0) 89 4623 17-0

**Italy:** +39 0 2 57603050

**UK:** +44 (0) 20 8204 8142

If an instrument or component must be returned, the method described on the following page should be followed to expedite servicing and reduce your downtime.

## Return Authorization

---

All instruments and components returned to the factory must be accompanied by a Return Authorization Number issued by our Customer Service Department.

To issue a Return Authorization number, we require:

- The model and serial number of the instrument
- A list of items and/or components to be returned
- A description of the problem, including operating settings
- The instrument user's name, mailing address, telephone, and fax numbers
- The shipping address for shipment of the instrument to you after service
- Your Purchase Order number and billing information for non-warranty services
- Our original Sales Order number, if known
- Your Customer Account number, if known
- Any special instructions